


LASER HAIR REGROWTH BRUSH

1. Comb Teeth
2. Lasers
3. RF / EMS
4. On / Off
5. Photon
6. Massage, Laser / Red Photon
7. Charging Indicator
8. Charging Slot
9. Base
10. Charging Cable

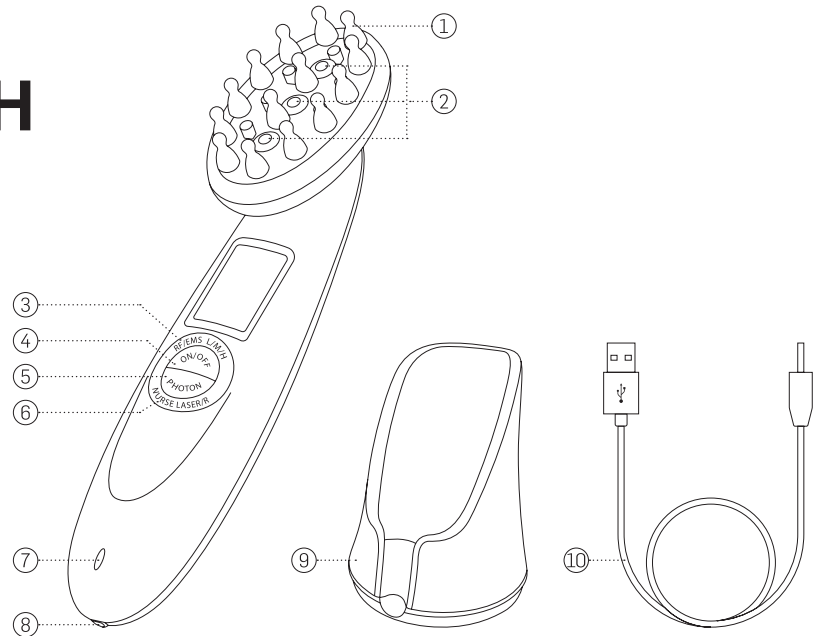

Laser Hair Regrowth Brush uses micro vibration, radiofrequency therapy, photon therapy, laser therapy, and electromagnetic stimulation therapy. These processes combined promote blood circulation and metabolism, and stimulates hair follicles which promotes hair growth and reverses hair loss. This device is suitable for those losing hair due to aging, hormones, or immune disorders.

How To Use:

1. Ensure the Laser Hair Regrowth Brush is charged and press the “On/Off” button.
2. Press the “Photon” button to adjust setting.
3. Gently run the brush through your hair. Be sure to make as much contact with your scalp as possible.
4. Use for no more than 10 minutes every day to ensure the best possible result.
5. To turn off, hold down the “On/Off” button.

Functions:

1. Micro-Vibration – Press the “Photon” button once to access this function.
2. Photon – Press the “Photon” button twice to access this function. There are three Photon settings, the device will cycle through them with each successive tap of the “Photon” button.
3. RF/EMS – The Electromagnetic Stimulation (EMS) and Radiofrequency (RF) functions can be accessed by pressing the “On/Off” button once while the device is on. To adjust the power level, tap the “On/Off” until the desired setting is reached.

Safety Notes:

- Do not use this device if you have any implanted metal or electrical medical instruments.
- Do not disassemble Laser Hair Regrowth Brush. This may cause damage to the device.
- Cease use if Laser Hair Regrowth Brush starts to emit any abnormal sounds, smells, or heat.
- Use Laser Hair Regrowth Brush only for its intended purpose.
- Do not use any longer than the recommended time.


Please note that this device is not suitable for those with the following conditions:

- Heart disease, artificial heart or lungs, installed pacemaker
- Any implanted metal medical device
- Cancer
- Hypotension, fever
- Sunburned
- Pregnant, lactation period
- Skin disease
- Undergoing medical treatment

Before and After Use with the Laser Hair Regrowth Brush:



Suitable for the Following Conditions:

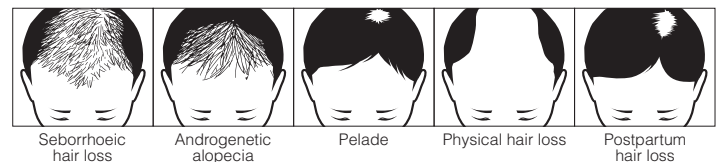

Specifications:

- Power Supply: AC100-240V, 50/50HZ
- Power: 3W
- Photon Wavelength: 630-650nm
- Size: 18.5*4.5*5.5cm
- Weight: 150g